

MiniTell

5 TIMES DAILY VIBRATING ALARM AND PILL BOX



OPERATING INSTRUCTIONS

Features:

- Five alarms per day either audio or vibration or both. Snooze function.
- Five compartments to hold daily medication (capacity - approx 10 ml.)
- Compact design for easy carriage in pocket or bag
- Locking switch to prevent accidental alteration of settings
- Large LCD display with LED backlight.

ISO9001:2000 and Medical CE (EUROPE) certified. RoHS approved

BATTERIES

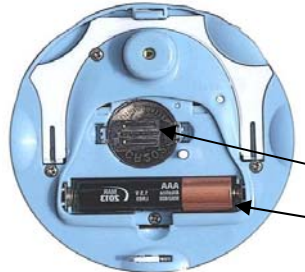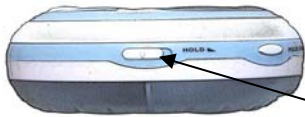

The MiniTell pill box is powered by one CR2025 lithium battery (preinstalled) for time function and one AAA battery for vibration and alarm sound function..

To install the AAA battery:


With the digital display facing upwards slide the '**OPEN**' switch on the rim to the right and separate the top body from the bottom.

Remove the battery case cover (picture shows cover already removed). Insert the AAA battery making sure it is placed with the correct polarity as indicated by the polarity symbols in the battery case recess. Replace battery case cover. Click bottom and top of unit back together.

Lithium CR2025 battery

AAA Battery

FILLING

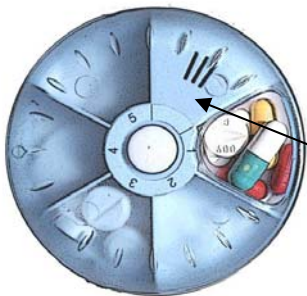

Turn the pill box so that the side with the pill compartments is facing you. Turn the translucent cover 1/6th round to obtain access to each compartment.

Fill the compartments numbered 1 to 5 with the medication corresponding to your prescription. i.e. Compartment 1 for first dose, compartment 2 for second...etc

Finish by turning the translucent cover so that all compartments are covered and check that all medication is in the correctly numbered compartment.

SETTING THE TIME



Turn the pill box so the that the digital display is facing you.

Slide the 'Hold' button (on the rim) to the left, to unlock settings.

Press the T-SET button for four seconds until **T-SET** appears on the LCD. Release button. The hours digits will flash.

Turn the Time Adjustment Ring (white serrated rim) either clockwise or anticlockwise to set the correct hour.

Press & release T-SET to confirm and the minutes digits will flash. Rotate the Time Adjustment Ring in the same way to set the correct minute.

Press & release T-SET again to confirm.

