



www.pivotell.co.uk



2026


PROGRAMMING ON A SMARTPHONE

Pivotell Smart

Important:

**Read the dispenser programming
instructions in full to understand
functionality first.**

HTTP mode and Wi-Fi Connection

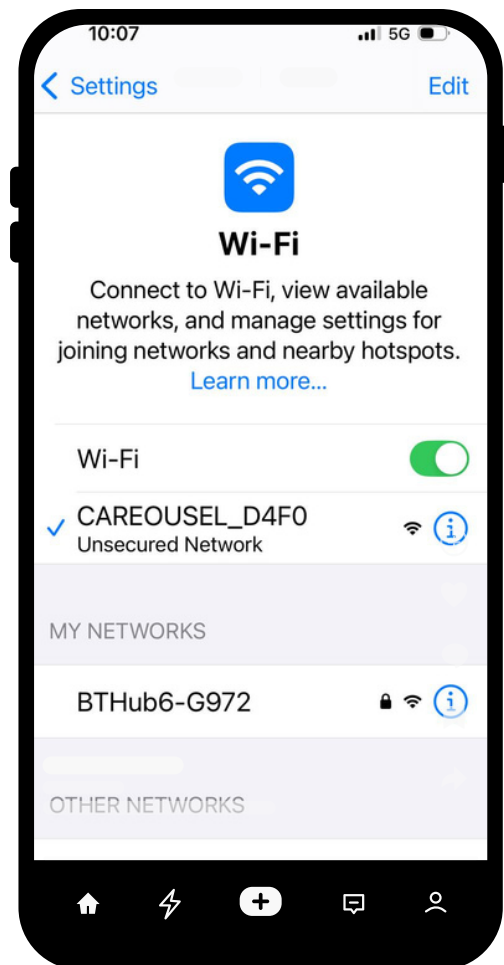

Step 1: On the Pill Dispenser

Press and hold buttons 1 and 3 together.
HTTP will display on LCD screen.

Step 2: Connect Smartphone to Dispenser Wi-Fi

- Open Wi-Fi settings on your smartphone
- Select CAREOUSEL_XXYY
- If not visible, turn Wi-Fi off and on, then try again



Access Settings

If you are reading this **on your smartphone**, tap here to open Settings:

Open Settings




Select **Settings** to program the device. Ignore the other options.

If you are reading this on another device or a printed copy:

Use your smartphone to:

- Scan the QR code, **or**
- Visit **192.168.4.1** in your phone's web browser

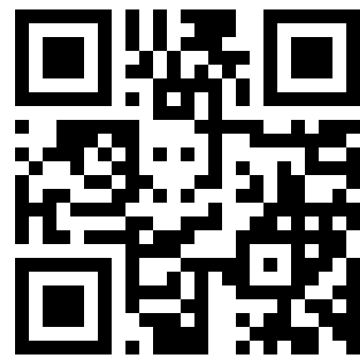


Refer to the programming manual for settings details. Press SAVE when finished. Setup is now complete. No further steps are required.



Tip: Share this instruction pdf and the weblinks to your phone Home Screen as it will make it much easier to return to it later.

Share icons look like this:

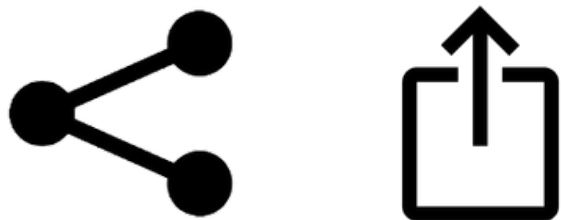

**Don't hesitate
to contact us for support**



office@pivotell.co.uk



www.pivotell.co.uk



01799 550979